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Instructions for use

Vacuum Activator silencos and silencos kids

Please read these instructions carefully before you use the product.

The manufacturer disclaims any liability for damages resulting from non-observance of these instructions.

1. Product description

The silencos vacuum activator consists of an elastic oral vestibule shield and an elastic membrane, a hard plastic lip attachment and attachable components (membrane funnel, stopper). For biofunctional early child (silencos kids) and adult (silencos) therapies. 1 piece!

2. Indication

- Interferences of the orofacial system's rest position
- Habitual mouth breathing
- Childlike habits (e.g. open mouth posture, lip biting, sucking)
- Primary snoring during habitual mouth breathing with velar genesis
- Dysgnathia due to habitual orofacial dysfunction
- Rehabilitation in terms of reestablishment of the self-contained rest position following tumor surgical measures and neurologic diseases involving orofacial functions
- Interferences of the oral phase of the swallowing process
- Xerostomia of various genesis
- Accompanying therapy in case of periodontal diseases to encourage salivation

3. Contraindication

Neurologic, pediatric and psychiatric disorders during which, a proper application is not possible.

4. Hazard and safety warnings

- Children should use silencos only while supervised by adults
- Consult your physician regarding night-time wear, as small parts may be swallowed or aspirated.
- See your treating physician immediately in case of periodontal problems (pressure spots, gum bleeding, pain).
- Ensure that individual parts are not interchanged to other patients to exclude cross-contamination.

5. Storage and shelf life information

Store in a dry and dust-free place.

6. Processing

The silencos vacuum activator consists of an elastic oral vestibule shield and an elastic membrane, a hard plastic lip attachment and attachable components (membrane funnel, stopper).

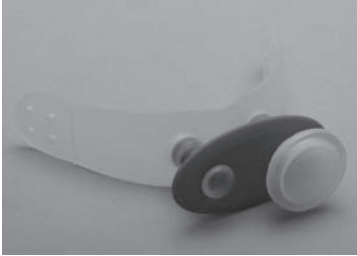
Fig. 1 shows the completely assembled vacuum activator

Fig. 2 shows the vacuum activator in its individual components (a) oral vestibule shield, b) lip attachment, c) membrane funnel, d) membrane, e) stopper).

The device has two openings at the lip attachment's plate that protrudes from the mouth: The membrane funnel is inserted into one opening; the other opening is closed with the stopper. This opening is used to attach, if necessary a medical tube system or pressure gauge system. The systems must be hermetically closed and pressure-tight coupled. The charged membrane is placed onto the funnel to create a pressurized area. Suction is used to check whether the membrane pulls inwards. If necessary, the (possibly damaged) membrane is replaced by a new one.

7. Application

The oral vestibule shield is inserted between the rows of teeth and the lips, whereby the lip attachment should be placed between the upper and the lower lip; the lips should be loosely closed. Prior to the initial insertion of the device examine the area between lips, cheeks, and tooth sectors for irregularities such as deep attached frenula or inconveniently shaped oral vestibules. In case the device compresses the frenulum or oral mucosa, the edges may be adequately trimmed with scissors. While training, the patient should not feel any pressure from the edge of the vacuum activator.

If an interdental tongue position is to be anticipated when activating silencos kids in early treatment of children, it is recommended to have a customized tongue blocker (see fig. 3) produced by the dental-technical laboratory. The tongue blocker is connected to the front openings of the mouth attachment piece by the dental staff and is positioned at the level of the occlusion plane in situ.

Figure 3.
Tongue blocker, customized

8. Therapeutic basic exercise

silencos offers various possibilities for training and therapy. It is mainly used as a training device during dental, orthodontic or logopedic treatments. The basic exercise of the biofunctional therapy with silencos consists of the initiation of the self-contained rest position of the orofacial systems. During that process, a negative pressure in the mouth creates a balance of forces between the outer soft tissues (lips and cheeks) and the tongue inside the mouth. At the same time, a tight closure between tongue and velum takes place. A negative pressure is created by saliva collection and the swallowing that follows. Once the saliva is swallowed, the patient breathes calmly through the nose and thus maintains the negative pressure. An active muscle tension is not necessary. The position of the silencos membrane gives information about the pressure conditions within the mouth. During this exercise, the membrane should be pulled into the tunnel to indicate the proper adoption of the self-contained rest position. Training should be conducted for 15 – 30 minutes several times a day. Aim of the basic exercise is to acquire and permanently maintain the self-contained state of rest during the phase the mouth is at rest. An excessive creation of negative pressure or numerous sole suction processes without swallowing should be avoided. Please consult your therapist regarding training frequency and duration, or further therapeutic measures with the silencos device.

9. Technical data

See material safety data sheet

10. Used symbols

- Lot number
- Please observe instructions for use

11. Additional information

- Insertion by specialized personnel only.
- Use only original attachment components. Please see your treating physician in case of loss.

12. Hygiene measures and storage

silencos and silencos kids is intended for the treatment of one patient only. Disassemble the device after use and remove visible dirt (e. g. with a soft toothbrush) and clean in lukewarm water with a detergent additive. Rinse with clear water and store the dry device in its storage box until the next use. Disinfect the device with an adequate disinfectant in case the patient suffers from a contagious disease. Additional information under www.silencos.de or www.silencos-kids.de.

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GmbH & Co.KG · Weissenhorner Str. 2 · 89250 Senden · Germany